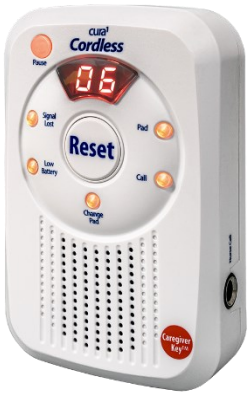
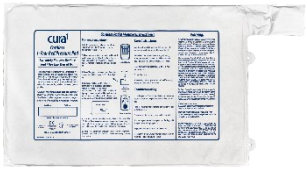
• Monitor up to 7 Pads or Mats simultaneously • Works with a single corded Premium or Universal Pad/ Mat, and up to six cordless Pads/Mats • Durable ABS plastic enclosure, and soft touch removable protection case and Caregiver Reset Key • Alarm memory and auto re-arm • Amber LED visual alarm indicator • Caregiver secure reset option • Battery operated or optional external plug pack • Up to 6 months battery life (nominal) • Easy to access battery compartment • Can be connected to a nurse call system • 3 selectable alarm tones • Easy thumbwheel adjustable alarm volume • "Pause" feature allows caregiver to attend the patient without setting the alarm off • Alerts for lost nurse call connection, lost AC power supply, lost transmitter signal & low battery		Self Reset	**2700 6Ch. Cordless Falls Monitor Compatability Matrix** Suitable for private residences and care homes, its versatility turns a basic room into a highly optioned care and monitoring facility offering multiple solutions in one simple to operate unit, our most powerful Falls Monitor yet. **What's in the box:** 1 x 2700 Cordless Monitor 1 x Quick Start Guide 1 x 2642 Caregiver Reset Key 1 x Durable Soft Case 1 x Non-Absorbent Mounting Strap 3 x High Quality AA Batteries	
				24 Month Warranty
3515	Cordless Under Mattress Bed Pad 1 Year (830 x 260)	✓		Place **UNDER THE MATTRESS** Sends **OUT OF BED** message when getting out of bed. Resets **OUT OF BED** when back in bed or via 2700 Falls Monitor ****Requires 2719 Transmitter to function****
2629	Cordless Bed Pad 1 Year (760 x 500)	✓		Place **ON TOP OF THE MATTRESS**. Sends **OUT OF BED** message when getting out of bed. Resets **OUT OF BED** when back in bed or via 2700 Falls Monitor. ****Requires 2719 Transmitter to function****
2630	Cordless Chair Pad 1 Year (380 x 250)	✓		Sends **OUT OF CHAIR** message when getting out of chair. Resets **OUT OF CHAIR** when back on chair or via 2700 Falls Monitor. ****Requires 2719 Transmitter to function****
2624	Cordless Floor Mat 1 Year (1220 x 620)			Position in an area of **FALLS RISK**. Initiates alarm when **stepped on**. Reset the alarm via 2700 Falls Monitor. ****Requires 2719 Transmitter to function****
2719	Cordless Pad Transmitter			Required by Cordless Pads / Mats. Sends a signal to the 2700 Falls Monitor. ****1 year battery life. Replace the transmitter after one year****

All radio frequency (RF) signal transmission and reception distances quoted are indicative only and based on testing carried out under ideal conditions. Actual transmission and reception distances may vary and be dependent on specific environmental conditions.

2690	Universal Single Fold Bed Pad 18 Months (760 x 255)	✓		Place **ON TOP OF THE MATTRESS** Sends **OUT OF BED** message when getting out of bed. Resets **OUT OF BED** when back in bed or via 2700 Falls Monitor. ****Requires 2720 Transmitter to function****
2689	Universal Under Mattress Bed Pad 18 Months (830 x 260)	✓		Place **UNDER THE MATTRESS** Sends **OUT OF BED** message when getting out of bed. Resets **OUT OF BED** when back in bed or via 2700 Falls Monitor. ****Requires 2720 Transmitter to function****
2691	Universal Chair Pad 18 Months (380 x 255)	✓		Sends **OUT OF CHAIR** message when getting out of chair. Resets **OUT OF CHAIR** when back on chair or via 2700 Falls Monitor. ****Requires 2720 Transmitter to function****
2692	Universal Floor Mat 18 Months (1220 x 610)			Position in an area of **FALLS RISK**. Initiates alarm when stepped on. Reset the alarm via 2700 Falls Monitor. ****Requires 2720 Transmitter to function****
2720	Universal Cordless Pad Transmitter			Required by Universal Pads / Mats if not connected via lead. Sends a signal to the 2700 Falls Monitor. ****18 month battery life. Replace the transmitter after 18 months****
2670	Cordless PreAlert Bed Pad 1 Year (830 x 260)	✓		Located on **TOP OF THE MATTRESS**. Sends **OUT OF BED** message when getting out of bed. Resets **OUT OF BED** when back in bed or via 2700 Falls Monitor. ****Transmitter Included****
2671	Universal Single Fold Floor Mat 18 Months (900 x 600)			Position in an area of **FALLS RISK**. Initiates alarm when stepped on. Reset on 2700 Falls Monitor. Can be folded for storage or transportation. ****Requires 2720 Transmitter to function for cordless operation****
2672	Universal Dual Fold Floor Mat 18 Months (1190 x 600)			Position in an area of **FALLS RISK**. Initiates alarm when stepped on. Reset on 2700 Falls Monitor. Can be folded for storage or transportation. ****Requires 2720 Transmitter to function for cordless operation****

All radio frequency (RF) signal transmission and reception distances quoted are indicative only and based on testing carried out under ideal conditions. Actual transmission and reception distances may vary and be dependent on specific environmental conditions.

3516	Premium Under Mattress Bed Pad 1 Year (830 x 260)	✓		Position **UNDER MATTRESS** Sends **OUT OF BED** message when getting out of bed. Resets **OUT OF BED** when back in bed or via 2700 Falls Monitor.
2626	Premium Bed Pad 1 Year (760 x 500)	✓		Position **ABOVE MATTRESS** Sends **OUT OF BED** message when getting out of bed. Resets **OUT OF BED** when back in bed or via 2700 Falls Monitor.
2628	Premium Chair Pad 1 Year (380 x 250)	✓		Position **ON CHAIR / WHEELCHAIR**. Sends **OUT OF CHAIR** message when getting out of chair. Resets **OUT OF CHAIR** when back on chair or via 2700 Falls Monitor.
2633	Premium Heavy-Duty Floor Mat 1 Year (1220 x 610)			Position in an area of **falls risk**. The 2700 Falls Monitor alarms when stood on. Reset via 2700 Falls Monitor.
2622	SafeTime Bed Pad 90 Days (760 x 500)	✓		Position **ABOVE MATTRESS** Sends **OUT OF BED** message when getting out of bed. Alarms when the pad **expires**. Resets **OUT OF BED** when back in bed or via 2700 Falls Monitor.
2621	SafeTime Chair Pad 1 Year (380 x 250)	✓		Position **ON CHAIR / WHEELCHAIR**. Sends **OUT OF CHAIR** message when getting out of chair. Alarms when the pad **expires**. Resets **OUT OF CHAIR** when back on chair or via 2700 Falls Monitor.
2721	iBeam			Sends alarm when triggered. Features a narrow detection zone to reduce false alarms, triggers when a patient breaks the beam. Reset on 2700 Falls Monitor.
2588	PET (Personal Emergency Transmitter)			Initiates alarm when **RED Button** is pressed. Reset on 2700 Falls Monitor.
2595	Wristband Alarm Transmitter			Staff or patients wear on the wrist. Use as a duress alarm or call button. Initiates alarm when **RED Button** is pressed. Reset on 2700 Falls Monitor.

All radio frequency (RF) signal transmission and reception distances quoted are indicative only and based on testing carried out under ideal conditions. Actual transmission and reception distances may vary and be dependent on specific environmental conditions.

3502	ActiveCare Wireless Handset (Includes wall bracket)	✓		Initiates alarm when **RED button** is pressed. Includes **GREEN** wireless reset button or reset via 2700 Falls Monitor.
2599	Wireless SOS Call Button			Wall mountable emergency call point. Initiates alarm when **RED Button** is pressed. Reset via 2700 Falls Monitor.
3503	Wireless Door Contact			Sends alarm when door is **opened**. Can also be used as a simple **push button** call point. Reset on 2700 Falls Monitor.
2698	Universal Transmitter	✓		Turn any third-party alarm device into a wireless transmitter. Can also be used as a wireless **door contact**. Togglable **Reset** button.
2558	Wireless Waterproof Call Button w/pullcord	✓		Initiates alarm when **RED button** is pressed or **cord is pulled**. Includes wireless **reset button** or reset on **2700 Falls Monitor.**
2406	cura1 6V Plug Pack			Optional **6 Volt DC** plug pack.
2642	Caregiver Reset Key			Allows caregivers to **reset alarms**. Allows authorised staff to **program** various cura1 equipment including **3500R Monitors** Key chain can be **attached** to belt or clothing. **Lightweight**.
3504	Mounting Bracket			2700 Falls Monitor slides into bracket. Bracket mounts onto surfaces with **double-sided tape** or **M3 screws**.
3501	ActiveCare LCD Pager	✓		Reports all alerts on a limited alpha / numeric 7-character screen. Can be mounted on a wall or carried by a caregiver. Fits on 3504 Bracket Mounts. ***For a soft cover see 3506****

All radio frequency (RF) signal transmission and reception distances quoted are indicative only and based on testing carried out under ideal conditions. Actual transmission and reception distances may vary and be dependent on specific environmental conditions.

2730	SnapOn Bed Pad 1 Year (800 x 255)			Place **ON TOP OF THE MATTRESS** Sends **OUT OF BED** message when getting out of bed. Resets **OUT OF BED** when back in bed or via connected monitor. ****Requires 2732 SnapOn Pad Transmitter or 2733 SnapOn Pad Cable to function****
2731	SnapOn Chair Pad 1 Year (420 x 180)			Sends **OUT OF CHAIR** message when getting out of bed. Resets **OUT OF CHAIR** when back in bed or via connected monitor. ****Requires 2732 SnapOn Pad Transmitter or 2733 SnapOn Pad Cable to function****
2732	SnapOn Pad Transmitter			Required by SnapOn Pads. Sends a signal to the 2700 Monitor. ****Replaceable CR2045 battery.****
2733	SnapOn Pad Cable 2m			Easily connects to the 2730 and 2731 SnapOn Pad range. Sends a signal to wired monitors through an RJ-10 connetor.
2669	Cordless Under Mattress PreAlert Bed Pad 1 Year (820 x 240)			Place **UNDER THE MATTRESS** Sends **OUT OF BED** message when at edge of bed. Resets **OUT OF BED** when back in bed or via 2700 Falls Monitor

All radio frequency (RF) signal transmission and reception distances quoted are indicative only and based on testing carried out under ideal conditions. Actual transmission and reception distances may vary and be dependent on specific environmental conditions.